

Disposable Laryngoscope for Anesthesia and Emergency Care

DUOSCOPE™ = 2 blades in 1

The next generation disposable emergency laryngoscope designed for the most popular sizes;


- Macintosh profile for adults
- Miller profile for pediatrics

High quality ABS plastic molding provides the strength and rigidity needed with greater patient safety than all-metal or hybrid blades.

- Gentle touch on the patient's teeth

The disposable DUOSCOPE is a cost effective premium intubation device for anesthesia, emergency situations and crash cart or kit.

- Eliminates cost and time associated with reprocessing blades
- Minimizes risk of cross-contamination during intubation

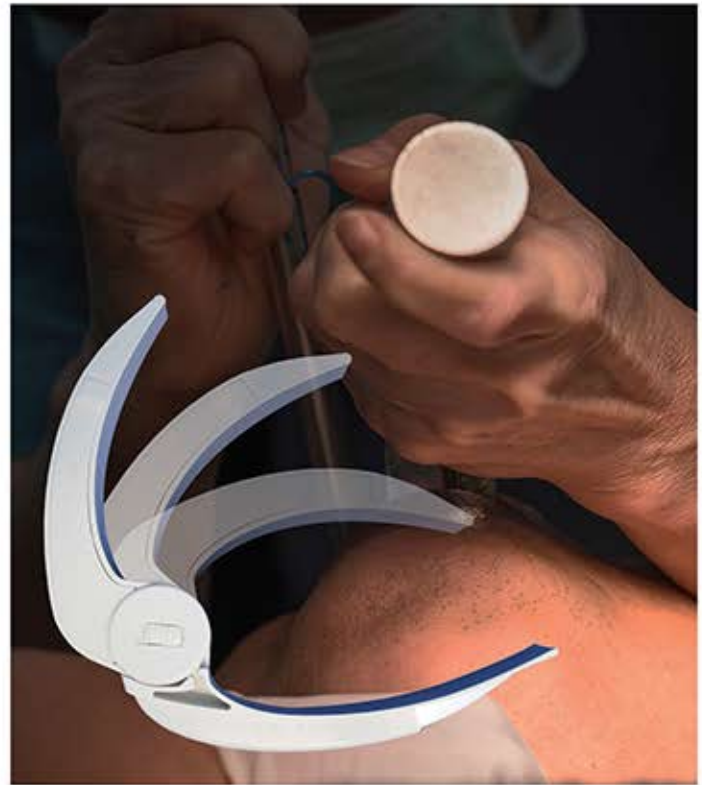

Trauma Rescue



Ambulance



Medical Air Service



Emergency Carts



Military



Anesthesiology

The DUOSCOPE's functional design facilitates immediate accessibility and instantaneous response in emergency situations.

No Assembly - Quick Action

Available in both Adult/Child and Adult/Adult dual blade configuration for convenience and ease of use.

Light-weight, "Two in One" blade without the need to carry multiple handles, blades or batteries.

Meets ISO standard for blade structural integrity, tested with an applied load of 40 kg force, there was no breakage for the full range of the blade.

Specifications

DUOSCOPE

Description

- Type - Disposable
- Material - ABS X-17 Plastic
- Blade Configuration / Size (closed position)

- Miller Size #1 + Macintosh Size #3	139mm H x 79mm W x 27mm D
- Miller Size #1 + Macintosh Size #4	169mm H x 79mm W x 28mm D
- Macintosh Size #3 + Macintosh Size #4	169mm H x 87mm W x 28mm D

- Dual Illumination Lamps, over 24 hrs. continuous illumination
- Light - LED, Super Bright Pure White Water Clear
- Aperture Opening - 90 degrees
- Battery - 3.0V Lithium-Manganese Coin Cell (2), easily removed for disposal/recycle
- Storage Conditions -22 degrees F to 140 degrees F (-30 degree C to 60 degree C)
- Non-sterile or sterile options
 - Each blade and handle packaged in a sealed pouch ready for ETO sterilization
- Meets International Standards for Laryngoscopes:
 - ISO: 7376:2009: Anesthesia and respiratory equipment - Laryngoscopes for tracheal intubation

3 Blade Configurations



Miller 1 +
Macintosh 3
RB2000



Miller 1 +
Macintosh 4
RB3000



Macintosh 3 +
Macintosh 4
RB4000



Airon

Airon Corporation

751 North Drive, Unit 6

Melbourne, Florida 32934

TEL +1.321.821.9433 / 888.448.1238 USA

FAX +1.321.821.9433



For more information visit
www.AironUSA.com